Saturn[®]

ST-PS0283

BATHROOM SCALES

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Bathroom Scales With Body Fat Analyzer

Dear Buyer!

We congratulate you on having bought the device under trade name "Saturn". We are sure that our devices will become faithful and reliable assistance in your housekeeping.

Avoid extreme temperature changes. Rapid temperature change (e.g. when the unit is moved from freezing temperature to a warm room) may cause condensation inside the unit and a malfunction when it is switched on. In this case leave the unit at room temperature for at least 1.5 hours before switching it on. If the unit has been in transit, leave it indoors for at least 1.5 hours before starting operation.

Contents

General safety instructions

Symbols in these instructions for use

Children

Special Safety Instructions For This Machine

Handling batteries

Inserting the Batteries

Location of Controls

Technical Specifications

Information on Measuring Body Fat

Body Fat

The Importance of Water for our Health

Guideline Figures

Use

Weighing Function

Saving Personal Data

Measurement and Analysis

Display Error Messages

Cleaning

General safety Instructions

Read the operating instructions carefully before putting the appliance into operation and keep the instructions including the warranty, the receipt and, if possible the box with the internal packing. If you give this device to other people, please also pass on the operating instruction.

- The device is designed exclusively for private and not for commercial use.
- Only use the device for the intended purpose.
- Do not use a damaged device.
- Keep away from heat, direct sunlight, damp (never immerse in liquids) and sharp edges.
- Do not try to repair the appliance on your own. Always contact an authorized

technician.

- Pay careful attention to the following "Special Safety Instructions".
- Service life – 2 years

Symbols in these Instructions for Use

Important information for your safety is specially marked. It is essential to comply with these instructions in order to avoid accidents and prevent damage to the machine:

WARNING:

This warns you of dangers to your health and indicates possible injury risks.

CAUTION:

This refers to possible hazards to the machine or other objects.

NOTE:

This highlights tips and information.

Children

•In order to ensure your children's safety, please keep all packaging (plastic bags, boxes, polystyrene etc.) out of their reach.

WARNING!

Caution! Do not allow small children to play with the foil as there is a **danger of suffocation!**

- Children should be supervised at all times in order to ensure that they do not play with device.

Special Safety Instructions For This Machine

WARNING

If plastic surface is moist or covered with water droplets, please wipe the surface dry as there is a danger of slipping.

CAUTION:

- Handle the balance with care, as it is a precise measuring instrument. Do not prop it or jump on it.
- Do not load the appliance with more than 180kg to avoid damaging the device.
- If the scales are carried or not going to be used, please store them only in a horizontal position.
- The body fat analysis scale is not suitable for the following persons:
 - Children under 10 years and adults over 85 years
 - Persons with symptoms of fever, oedema, or osteoporosis

- Persons undergoing treatment for dialysis
- Persons who use a pacemaker
- Persons who take cardiovascular medication
- Pregnant women
- Sportspersons who perform over 5 hours of intensive sporting activities each day
- Persona with a heartbeat lower than 60 beats/minute

Handing batteries

WARNING:

- Do not expose the battery to high temperatures or direct sunlight. Never throw batteries into the fire. There is a danger of explosion!
- Keep batteries away from children. They are not a toy!
- Do not open the batteries by force.

WARNING:

- Avoid contact with metallic objects. (Rings, nails, screws et cetera) there is a danger of short-circuiting!
 - As a result of a short-circuit batteries may heat up considerably or even catch fire, this may result in burns.
 - For your safety the battery poles should be covered with adhesive strips during transport.
 - If a battery leaks do not rub the liquid into contact with the liquid they should be washed thoroughly and eyes rinsed out with clear water, if symptoms persist please consult a doctor.
4. Please turn socks inside-out when washing.

Inserting the batteries

NOTE:

- Where required, the supplied batteries have been protected for transport with plastic foil.
1. Remove this foil before inserting the batteries into the appliance.
 2. Open the battery compartment on the back.
 3. Insert 2 batteries of type AAA/R03. Please ensure that the polarity is correct (see the floor of the battery compartment)!
 4. Close the battery compartment again.
 5. If the device is not used for a prolonged period of time, please take the batteries out in order to prevent leakage of battery acid.
 6. If the message "LO" appears on the display, replace the batteries.

CAUTION:

- Different battery types or new and used batteries must not be used

together.

- Batteries are not to be disposed of together with domestic wastes. Dispose of the used batteries at the relevant collection bin or at your dealer's.

Location of controls

1. Contact surfaces
2. Display
3. DOWN-button ▼
4. UP-button ▲
5. SET-button

Information on measuring Body Fat

Display symbol	Meaning
Kg/lb/st	Body weight
P0-9	User memory space
	Male
	Female
CM(FT;In)	Body height
Age	Age
FAT	Body fat display
TBW	Muscle mass display
MUS	Muscle mass display
BONE	Bone weight display
%	Percentage-body fat/body water/muscle mass/bone weight

WARNING:

Do not use the scales if you use a pacemaker or another medical implant with electronic components. The function of the implant may be impaired the measurement.

NOTE:

The device can only function correctly if you stand on the contact surface with bare, dry feet!

The main principle of this device is based on the measurement of the electric impedance within the human body. In order to determine the actual

body date, the device reads a weak, unnoticeable electric signal through your body.

This signal is generally safe! However, please read our special safety instructions!

The method of measurement is referred to as "bioelectric impedance analysis" (BIA), a factor which is based on the relationship between body weight and associated body water percentages as well as on other biological data (age, gender, body weight).

There are heavy persons with a well-balanced percentage between fat and muscle mass, and there are seemingly slim persons with significantly high body fat values, Fat is a factor that puts serious strain on the body, in particular the cardio-vascular system. Therefore, controlling both values (fat and weight) is a prerequisite for health and fitness.

Body fat

Body fat is an important component of the body. From protecting the joints as a vitamin repository to regulating the body temperature, fat plays a vital role in the body. Thus the goal is not radically reduce the proportion of body fat, but rather to achieve and maintain a healthy balance between muscle and fat.

Guideline values for body fat and body water content as a percentage (approximate figures)

Age	Women		Men		
	Fat	Water	Fat	Water	
≤30	4.0-20.5	66.0-54.7	4.0-15.5	66.0-58.1	Underweight
	20.6-25.0	54.6-51.6	15.6-20.0	58.0-55.0	Normal
	25.1-30.5	51.5-47.8	20.1-24.5	54.9-51.9	Overweight
	30.6-50.0	47.7-34.4	24.6-50.0	51.8-34.4	Obese
	4.0-25.0	66.0-51.6	4.0-19.5	66.0-55.4	Underweight

>30	25.1-30.0	51.5-48.1	19.6-24.0	55.3-52.3	Normal
	30.1-35.0	48.0-44.7	24.1-26.5	52.2-49.2	Overweight
	31.1-50.0	44.6-34.4	26.6-50.0	49.1-34.4	Obese

Underweight mark: -
Normal Mark: 0
Overweight Mark: +
Obese mark: ++


Use

The device has a very sensitive sensor. If you want to adjust the settings, simply touch the buttons gently with your finger. If the device does not respond to touch, clean the glass surface with a soft, dry cloth.

NOTE:

The scales are factory set with the units "kg" and "cm". To change the settings to "lb" or "kg" and "Ft: in", proceed as follows: On the back side of the scale beside the battery compartment a switch can be found, push the below switch, choose between the measuring units kg/lb/st.

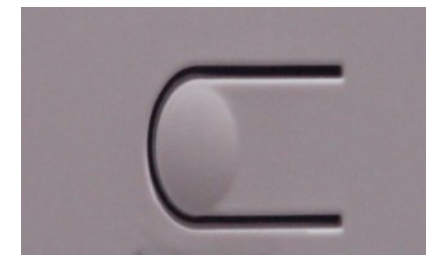Weighing Function

- Place the scale on a hard, flat surface (avoid carpets). An uneven floor contributes to inaccurate Measurements.
- Stand on the scales with both feet. The device switches on automatically.
- Distribute your weight equally and wait for the scale to determine your weight.

- The digits on the display will flash before your exact weight will be shown.
- In order to extend the battery life, the device will automatically turn off 10 seconds after you have stepped off the scale.

Saving Personal Data

- i NOTE:**
- The scale has default values that you can change, as described below.
 - Carry out the following steps in sequence and within a few seconds. Otherwise, the device will automatically cancel input mode.
 - Press the button SET. The memory space will flash on the display (USER0. Confirm the setting by pressing SET or select another memory space by pressing the UP/DOWN button. Then press SET to go to the next setting.
 - The symbol for "male" or "female" will flash. Confirm the setting by pressing SET or use the UP/DOWN button to choose the symbol for your gender. Then press to go the next setting.
 - The height setting will flash. Select your height by pressing UP or DOWN button. Holding down one of the buttons will accelerate the process. Confirm your entry by pressing SET.
 - The setting AGE will flash. Select your age by pressing UP or DOWN button. Holding down one of the buttons will accelerate the process. Confirm your entry by pressing SET.

- i NOTE:**
- The device can only function correctly if you stand on the device with bare, dry feet!
 - If necessary wipe the contact surface improve results.
 - Due to the significant fluctuation of body fat, it is different to attain a standard value. On average the optimal testing time is between 6 p.m. and 8 p.m., since the body fat remains very consistent during this period of time.
 - Press the SET button. The preset stored position date appears on the display (USER flash). Use the UP/DOWN button to select the stored position at which your date is saved.
 - Wait for the USER display and 0.0kg (0.0lb/st) to stop flashing.
 - Stand on the contact surfaces (1) of the scales with both feet.
 - Firstly your weight will flash on the display.
 - After a few seconds, the following results will be shown:

1. -the body fat content in % (FAT symbol)
-the body water content in % (TBW symbol)
-the results analysis in the form of a bar
- 2 -the muscle mass content in % (MUS symbol)
-the bone mass content in % (BON symbol)
-the calorie requirement per day
- 3 -The present gender
-The present type
-The present height
-The present age (AGE)

Display error Message

- i NOTE;**
- If an error message appears on the display, wait until the device is turned off. Turn it on again to restart the operation. If the error display has not disappeared, take the batteries out of the battery compartment for 10 seconds.

Display symbol	Meaning
8888	Self-test (when switching on) to check that the numeric display functions properly
----	The scales stand faultily
Lo	Battery is low
Err	The percentage of body fat content lies above or below the limits of the measuring scale
EEEE	The weight lies above the allowed maximum weight

Cleaning

- Clean the device with a slightly moist cloth. Make sure that no humidity or waters seeps into the device.
- Use a commercially sold detergent or soap for cleaning.
- Dry the device well after cleaning with a dry and soft cloth.

Technical Specifications

Model: ST-PS0283
 Power supply: batteries 1.5V
 Type "AAA"/R03"
 Net weight of device: 2KG
 Memory capacity: 10
 Height alignment: from 75 to 225 cm (2'5"-7'4 1/2 "Ft; in)
 Age: from 10 to 85 years
 Minimum weight: from 2.5kg (5.5lb/0; 4st)
 Maximum weight: 150Kg (330.7lb/23; 6st)
 Weight measurement unit: 100g increments
 Body fat measuring unit: 0.1%
 Body fat measuring range: 4.0%-50.0%
 Body water measuring range: 35.0%-85.0%
 Muscle mass measuring: 1.0%-50.0%

Subject to technical changes without prior notice!

This device has been tested according to all relevant current CE guidelines, such as electromagnetic compatibility and low voltage directives, and has been constructed in accordance with the latest regulations.

SET

Bathroom Scales	1
Instruction Manual with Warranty Book	1
Package	1

ENVIRONMENT FRIENDLY DISPOSAL

You can help protect the environment!
 Please remember to respect the local regulations: hand in the non-working electrical equipment to an appropriate waste disposal center.

The manufacturer reserves the right to change the specification and design of goods.

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GB

INTERNATIONAL MANUFACTURER'S WARRANTY

The Warranty is provided for the period of 2 years or more if the Law on User Right of the country where the appliance was bought provides for a greater minimum warranty period.

All terms of warranty comply with to the Law on Protection of the Consumer Rights and are governed by the legislation of the country where the product was bought. The warranty and free of charge repair are provided in any country where the product is delivered to by Saturn Home Appliances or authorized representative thereof, and where no restriction of import or other legal provisions interfere with rendering of warranty services and free of charge repair.

Cases uncovered by the warranty and free of charge repair:

1. Guarantee coupon is filled in improperly.
2. The product has become inoperable because of the Buyer's nonobservance of the service regulations indicated in the instruction.
3. The product was used for professional, commercial or industrial purposes (except for the models, specially intended for this purpose that is indicated in the instruction).
4. The product has external mechanical damage or damage caused by penetration of the liquid inside, dust, insects and other foreign objects.
5. The product has damages caused by non-observance of the rules of power supply from batteries, mains or accumulators.
6. The product was exposed to opening, repair or modification of design by the persons, unauthorized to repair; independent cleaning of internal mechanisms etc. was made.
7. The product has natural wear of parts with limited service life, expendables etc.
8. The product has deposits of scale inside or outside of fire-bar elements, irrespective of used water quality.
9. The product has damages caused by effect of high (low) temperatures or fire on non-heat resistant parts of the product.
10. The product has damages of accessories and nozzles which are included into the complete set of product shipment (filters, grids, bags, flasks, cups, covers, knives, sealing rings, graters, disks, plates, tubes, hoses, brushes, and also power cords, headphone cords etc.).

On the warranty service and repair, please contact the authorized service centers of "Saturn Home Appliances".

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GB

WARRANTY COUPON

**When purchasing the product,
please require its checking before you;
BE SURE
that the goods sold to you, are functional and complete
and that the warranty coupon is filled in correctly.**

This warranty coupon confirms the absence of any defects in the product you purchased and provides for free of charge repair of the product failed through the fault of the manufacturer throughout the period of warranty service and free of charge repair.

Failing the presentation of this coupon, in case of its improper filling in, infringement of factory seals (if any), and also in cases indicated in the warranty obligations, the claims are not accepted, and no warranty service and free of charge repair is made!
The warranty coupon is valid only in the original copy with the stamp of trading organization, signature of the seller, date of sale and signature of the buyer.

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